


Eleganza 3 LE

Medical Surgical Hospital Bed

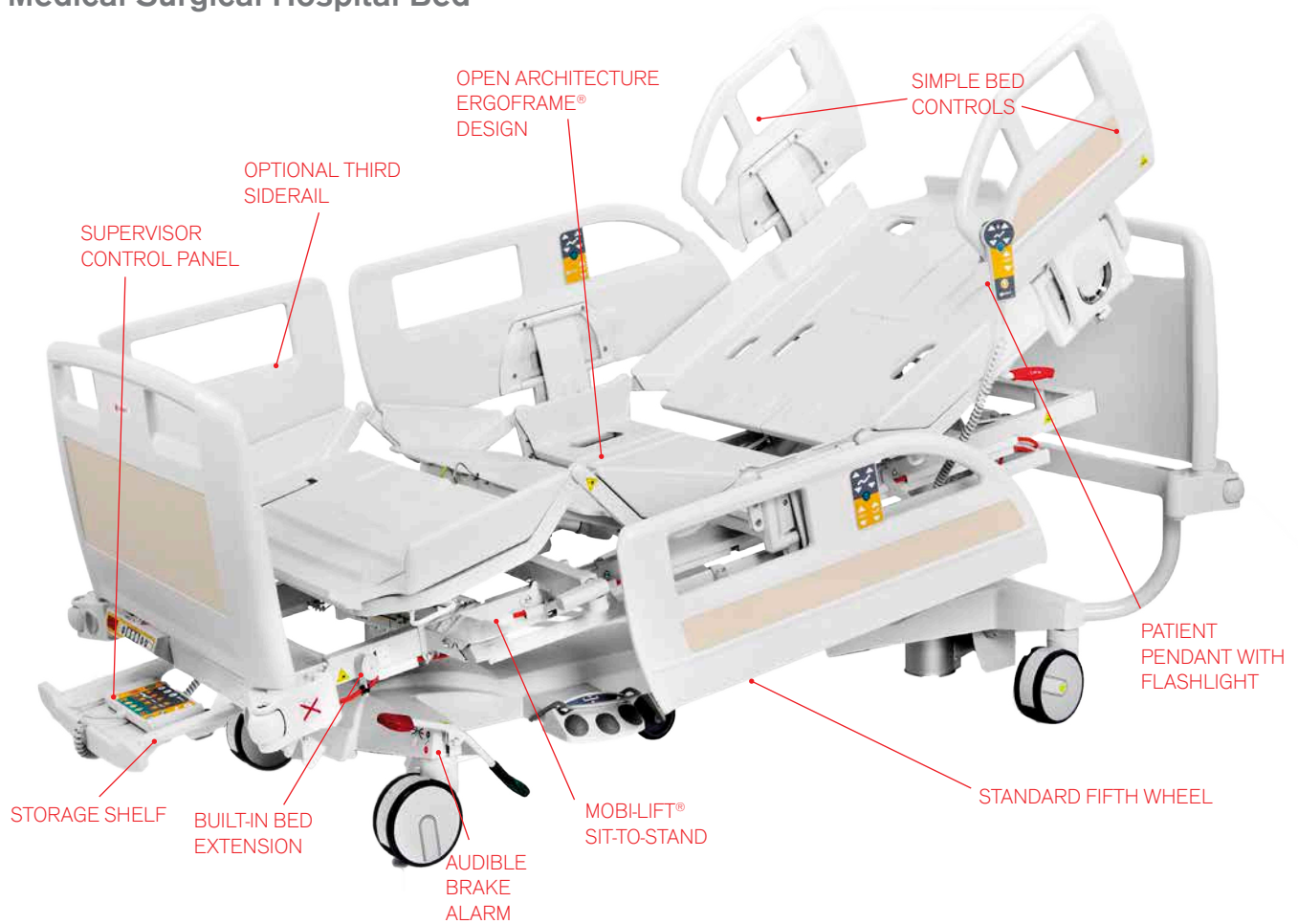

SPECIFICATIONS

Safe working load	550 lbs.
Frame height low	15.5 in.
Frame height high	30.5 in.
Frame width	41 in.
Frame length extended	97.7 in.
Frame length retracted	84.4 in.
Frame weight	315 lbs.
Trendelenburg	+/- 16 deg.
Max knee angle	30 deg.
Max HOB angle	70 deg.

STANDARD FEATURES

- 5-year warranty
- Mobi-Lift® sit-to-stand
- Integrated scale
- Multi-zone bed exit alarm
- Brake not set alarm
- Standard 5th wheel
- One-touch chair position
- One-touch bed egress position
- One-touch examination position
- One-touch auto-contour position
- Built-in bed extension
- Quiet DC motors with battery back-up
- 6 in. Tente® casters
- Removable sleep deck for easy cleaning and maintenance without tools

SURFACE OPTIONS

Eleganza 3 LE is designed with open architecture to allow flexibility to choose a surface that meets your specific wound care needs. Regardless of the surface, our Ergoframe® positioning system moves in harmony with the surface to help reduce pressure in the sacral area as the bed articulates.

Eleganza 3 LE



Mobi-Lift® Sit-to-Stand

- Built-in sit-to-stand device
- Handles on both sides of the bed
- Simple controls allow raising and lowering of bed
- Allows active participation from patient during mobilization



Integrated LED light

- Bright LED light for night time exams
- Integrated into patient pendant



5th Wheel Steering

- Enhances maneuverability of bed during transport



ErgoFrame®

- Helps reduce patient slide to foot of bed
- Auto-contour button raises head and knee simultaneously to maintain position and increase comfort



Open Architecture

- Choose a surface to meet your specific needs
- Removable deck for easy cleaning and maintenance
- Lightweight polymer material will not rust



Modular Design

- Flat deck, modular design is easy to clean and maintain
- Joint-free technology for improved infection control
- Top-down access for easy bed maintenance



Simple Controls

- Intuitive interface
- One touch chair
- One touch egress
- One touch auto-contour
- One touch examination position
- One touch CPR
- Lock-out functionality
- Safety GO button

Consider the Difference

With increased safety to prevent accidental bed movement, press any green "GO" button first to wake up the bed and adjust its position.

